

Retinal Image Selection Guide

แนวทางการเลือกภาพจอประสาทตาสำหรับแพทย์

1 Scope ขอบเขต

This guide sets out the properties a retinal image should have before you select it from hospital storage for review in Retinal Review Workbench. Grading and annotation are done later, inside the Workbench. Nothing needs to be graded, scored, or recorded at selection.


2 Image specification คุณสมบัติของภาพ

CRITERION	REQUIREMENT	NOT ACCEPTABLE
Modality ชนิดภาพ	Colour fundus photograph of the retina.	OCT, external-eye, anterior-segment, or other non-retinal images.
Source class การจัดกลุ่มต้นทาง	Already classified Refer or Not-refer in the hospital record.	Images with no existing classification.
Field of view ขอบเขตภาพ	Optic disc and macula both inside the image, with the retina filling most of the field.	Disc or macula cut off, or only a small part of the retina visible.
Focus ความคมชัด	Major vessels sharp. Small vessels around the disc and macula visible.	Blur that hides vessel detail.
Exposure ความสว่าง	Retina evenly visible, without large dark or washed-out areas.	Under- or over-exposure that hides the macula or a large area.
Artefacts สิ่งรบกวน	No major obstruction over the retina.	Eyelid, eyelash, reflection, or media opacity covering much of the field.
File ไฟล์	The original file as stored: JPEG, PNG, TIFF, or single-frame ophthalmic DICOM.	Screenshots, crops, annotated or edited copies, multi-frame or non-ophthalmic DICOM.
Eye (laterality) ตาข้างไหน	Left or right, if known from the record.	Guessing. If unsure, use Unknown .
Duplicates ภาพซ้ำ	Each image is selected once.	An image you know you have already selected.

3 Decision การตัดสินใจ

● Select ใช้ The image meets every criterion in Section 2.	● Review first ตรวจสอบซ้ำ The image is borderline on focus, exposure, or field. Select it only if you would still read it clinically.	● Do not select ไม่ใช้ Any item in the Not acceptable column applies, or the file cannot be opened.
--	---	---

4 Reference field and notes ภาพอ้างอิงและข้อควรทราบ

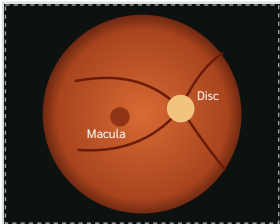

Adequate field (schematic drawing, not a patient image): disc and macula both inside the frame.

1. **Do not guess missing details.** Leave laterality or other unknown information as **Unknown**.
2. **Do not alter the file.** No renaming, cropping, annotating, converting, or re-saving before staging.
3. **Use only the hospital-approved staging process** to send selected images to the Workbench.
4. **Do not add patient identifiers** to file names or images.
5. **If more than one image of the same eye is available**, select the clearest one.