


Retinal Review Workbench - Clinician User Manual

This manual describes the final clinician interface. Screenshots are generated from the built application with public/synthetic data by `scripts/manual/capture_manual.py`.

1. Open the application

Open the review workstation at `http://127.0.0.1:8000/app/` after the operator starts `START.cmd`.

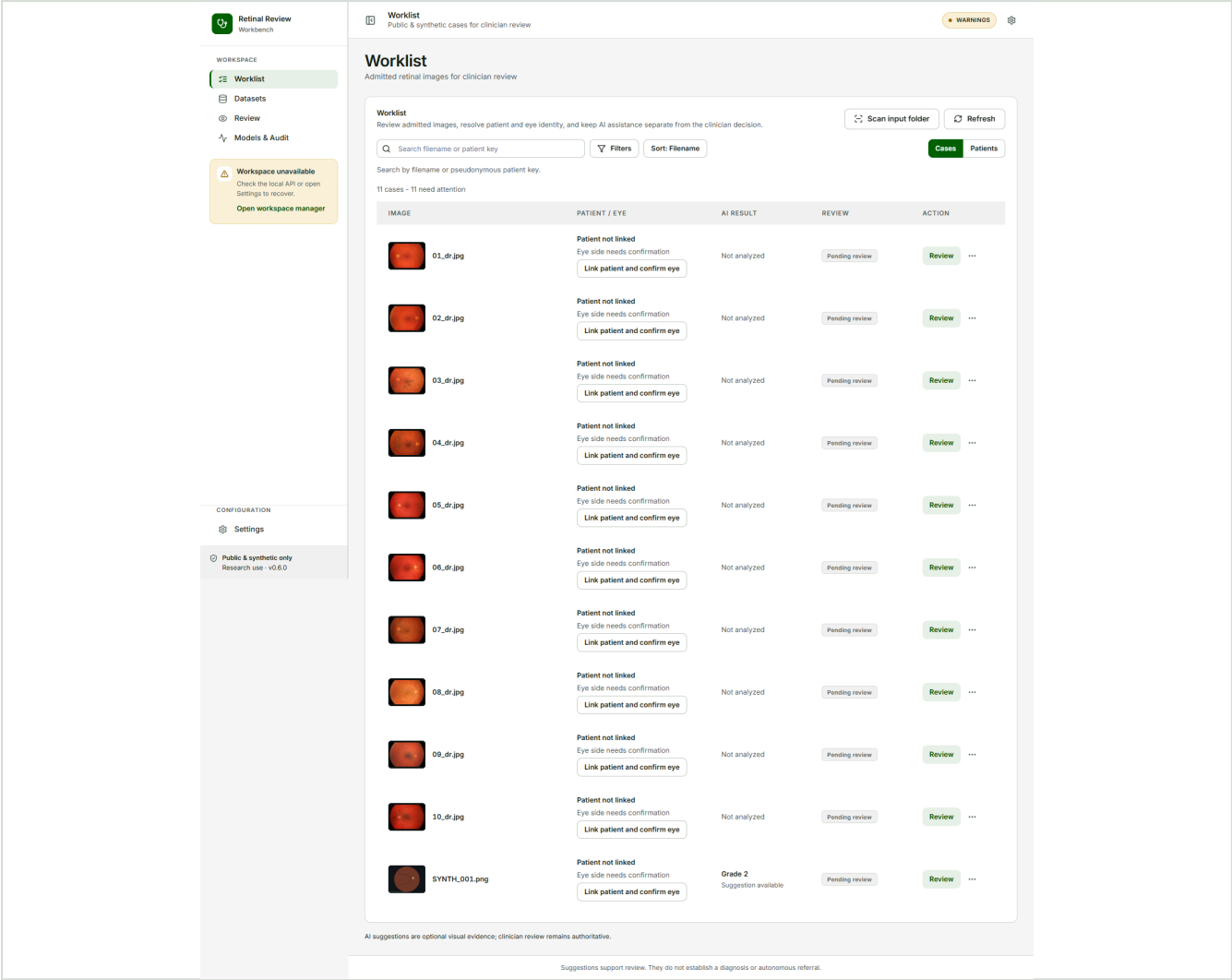

Worklist

2. Select a Workspace

Use Settings to create or select a Workspace. A Workspace defines the input folder, output folder, and SQLite review database. The application does not delete source images when a saved profile is removed.

Retinal Review
Workbench

WORKSPACE

- Worklist
- Datasets
- Review
- Models & Audit

Workspace unavailable

Check the local API or open Settings to recover.

Open workspace manager

CONFIGURATION

- Settings

Public & synthetic only

Research use - v0.6.0

Settings
Workspace, runtime, and clinician preferences

WARNINGS

Settings

Manage the local workspace used by this review station

Saved workspaces

Profiles keep local folders and the SQLite review database together.

+ New workspace

No workspaces saved

Create a profile to choose the local input folder, output folder, and SQLite database for this review station.

+ Create the first workspace

Workspace selection changes local review context only. It does not imply that files were moved, copied, uploaded, or sent to a model service.

AI Model Connection

Optional provider-neutral connection for a hospital LAN GPU server or local GPU API.

Not configured

Connection name *

Model API URL *

Model API

https://model-api.example

Access token (optional)

Enter only when required

Stored only by the running backend connection. It is never returned to this page or written to Workspace storage.

Test connection

Save

Suggestions support review. They do not establish a diagnosis or autonomous referral.

Workspace settings

3. Scan the input folder

Choose **Scan input folder** from Worklist. JPEG, PNG, TIFF, and supported ophthalmic DICOM can be admitted for review. MRI, OCT, PACS objects, and other non-fundus DICOM appear as **Unsupported modality**; they are not sent to the retinal models.

Retinal Review
Workbench

WORKSPACE

Worklist

Datasets

Review

Models & Audit

Workspace unavailable

Check the local API or open
Settings to recover.

Open workspace manager

CONFIGURATION

Settings

Public & synthetic only

Research use - v0.0.0

Worklist

Public & synthetic cases for clinician review

WARNINGS

Worklist

Review admitted images, resolve patient and eye identity, and keep AI assistance separate from the clinician decision.

Search filename or patient key

Filters

Sort: Filename

Scan input folder

Refresh

Cases

Patients

Search by filename or pseudonymous patient key.

11 cases - 11 need attention

IMAGE	PATIENT / EYE	AI RESULT	REVIEW	ACTION
01_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
02_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
03_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
04_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
05_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
06_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
07_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
08_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
09_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
10_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
SYNTH_001.png	Patient not linked Eye side needs confirmation Link patient and confirm eye	Grade 2 Suggestion available	Pending review	Review ...

AI suggestions are optional visual evidence; clinician review remains authoritative.

Suggestions support review. They do not establish a diagnosis or autonomous referral.

Input and worklist

4. Resolve patient and eye context

Use the Worklist identity action to confirm or set the pseudonymous patient key and laterality. Patient identity and eye laterality are independent fields. Leave a case unlinked when evidence is insufficient.

Retinal Review
Workbench

WORKSPACE

Worklist

Datasets

Review

Models & Audit

Workspace unavailable

Check the local API or open
Settings to recover.

Open workspace manager

CONFIGURATION

Settings

Public & synthetic only

Research use - v0.0.0

Worklist

Public & synthetic cases for clinician review

WARNINGS

Worklist

Review admitted retinal images for clinician review

Review admitted images, resolve patient and eye identity, and keep AI assistance separate from the clinician decision.

Search filename or patient key

Filters

Sort: Filename

Scan input folder

Refresh

Cases

Patients

Search by filename or pseudonymous patient key.

11 cases - 11 need attention

IMAGE	PATIENT / EYE	AI RESULT	REVIEW	ACTION
 01_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 02_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 03_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 04_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 05_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 06_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 07_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 08_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 09_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 10_dr.jpg	Patient not linked Eye side needs confirmation Link patient and confirm eye	Not analyzed	Pending review	Review ...
 SYNTH_001.png	Patient not linked Eye side needs confirmation Link patient and confirm eye	Grade 2 Suggestion available	Pending review	Review ...


AI suggestions are optional visual evidence; clinician review remains authoritative.

Suggestions support review. They do not establish a diagnosis or autonomous referral.

Patient and eye context

5. Open Review

Open an eligible case from Worklist. The retinal image remains the primary surface. Use Fit, 100%, zoom, pan, full-screen, and the original-pixel coordinate inspector as needed.

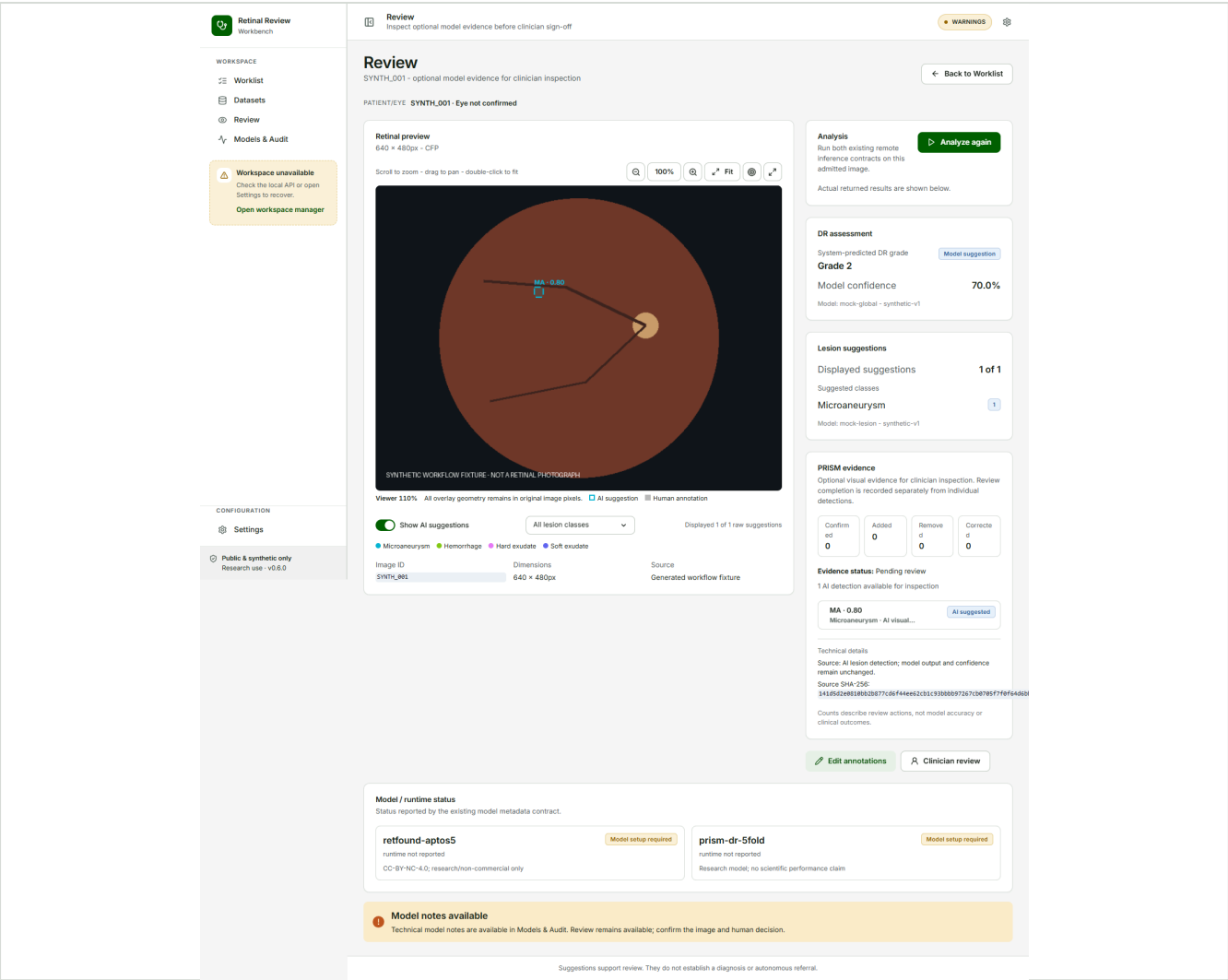


Review

6. Use optional model evidence

Select **Analyze** when model assistance is needed. RETFound provides a DR grade suggestion. PRISM-DR provides optional lesion ROI overlays. Toggle **Show AI suggestions** and use the lesion-class filter to reduce visual density. Each overlay displays its class and model score, for example `HE - 0.90`.

AI ROIs are evidence, not unresolved tasks. You may select an ROI to inspect it, edit or remove it when a recorded correction is needed, or leave it unchanged. You do not need to confirm or reject every detection before saving the clinician review.

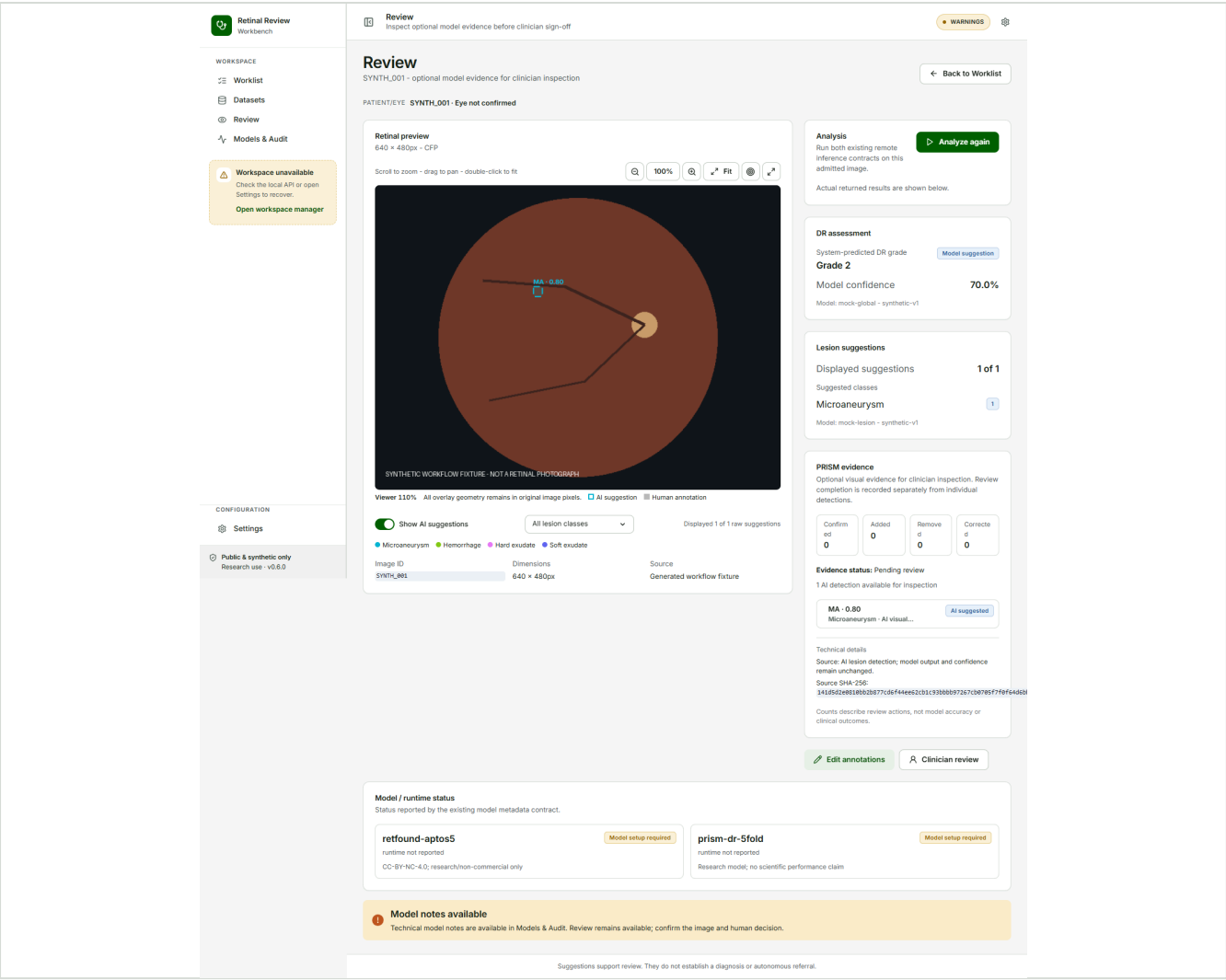

Optional overlays

7. Record clinician review

Open **Clinician review**, enter the reviewer name, and save an accepted, corrected, incorrect, or escalated decision. The final human review is authoritative and is independent of individual PRISM ROI actions.

Retinal Review
 Workbench

WORKSPACE

- Worklist
- Datasets
- Review
- Models & Audit

Workspace unavailable
 Check the local API or open Settings to recover.
[Open workspace manager](#)

CONFIGURATION

- Settings

Public & synthetic only
 Research use - v0.6.0

Clinician Review
 Human sign-off and review record

WARNINGS

Clinician Review
 SYNTH_001 - human sign-off and review record

Review summary
 AI output is a suggestion; this page records the clinician decision.

Scroll to zoom - drag to pan - double-click to fit

Viewer 102% All overlay geometry remains in original image pixels.

☒ AI suggestion
 ☐ Human annotation

AI-predicted grade	Human annotations	AI lesion suggestions
Grade 2	0	1

Current review state

Pending review

Clinician sign-off
 Save a human review action with an optional remark.

Reviewer name *

Final DR grade
 Select grade 0-4

Remark

Accept AI grade

Correct / set grade

Mark AI incorrect

Escalate

Review record
 Persisted in the existing SQLite case store.
 No clinician review saved yet.

Back to Review


Back to Worklist

Suggestions support review. They do not establish a diagnosis or autonomous referral.

Clinician review

8. Use dense annotation only when needed

Open ****Edit annotations**** for geometry-level work. Human annotations are separate from AI evidence and remain in original-image coordinates. CVAT Online is optional; use it only when a dense correction workflow is actually required.

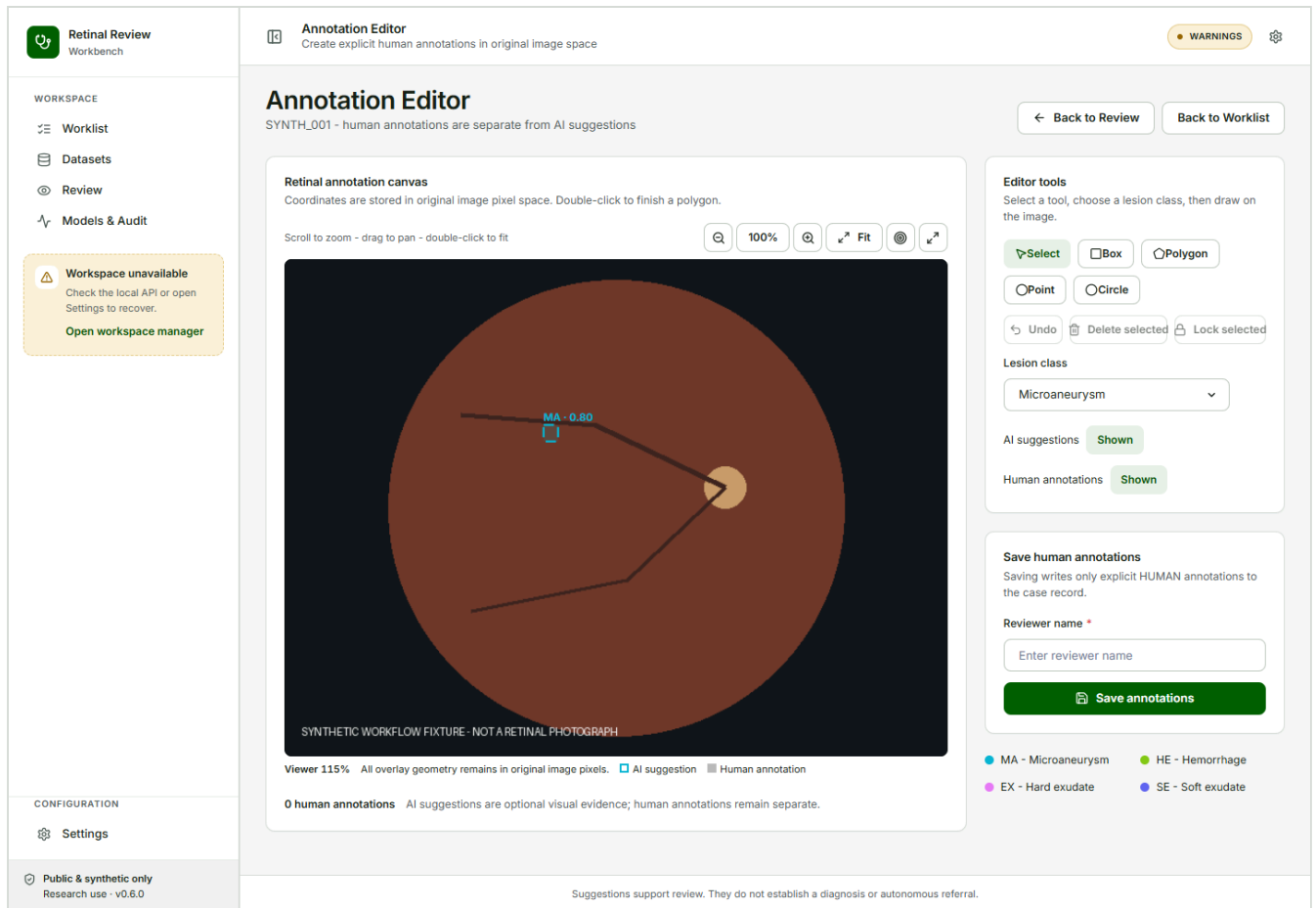

Annotation editor

9. Inspect Dataset and Models & Audit

Datasets shows the active Workspace, eligibility, annotation status, and export action. Models & Audit is read-only: it shows model IDs and revisions, inference metadata, detection counts by class, source and analysis SHA-256 values, and annotation provenance.

Retinal Review
Workbench

WORKSPACE

Worklist

Datasets

Review

Models & Audit

Workspace unavailable

Check the local API or open Settings to recover.

Open workspace manager

CONFIGURATION

Settings

Public & synthetic only

Research use - v0.6.0

Datasets

Manage retinal imaging workspaces

WARNINGS

Datasets

A reproducible manifest of the active Workspace, with AI and clinician provenance kept separate.

Export manifest

Dataset manifest

Export metadata only; source images and case records are never changed.

Refresh

Scope

No active Workspace

Images 11

Annotations 1

Training-ready 0

Open an active Workspace to export a dataset manifest.

All 11

Training-ready 0


Needs review 11

Excluded 0

IMAGE	PATIENT / EYE	CLINICIAN GRADE	ANNOTATIONS	DATASET STATUS
01_dr.jpg 01_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
02_dr.jpg 02_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
03_dr.jpg 03_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
04_dr.jpg 04_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
05_dr.jpg 05_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
06_dr.jpg 06_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
07_dr.jpg 07_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
08_dr.jpg 08_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
09_dr.jpg 09_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
10_dr.jpg 10_dr	Unlinked / Eye not confirmed	No final grade	None	Needs review
SYNTH_001.png SYNTH_001	Unlinked / Eye not confirmed	No final grade	1 AI	AI only

Suggestions support review. They do not establish a diagnosis or autonomous referral.

Dataset export

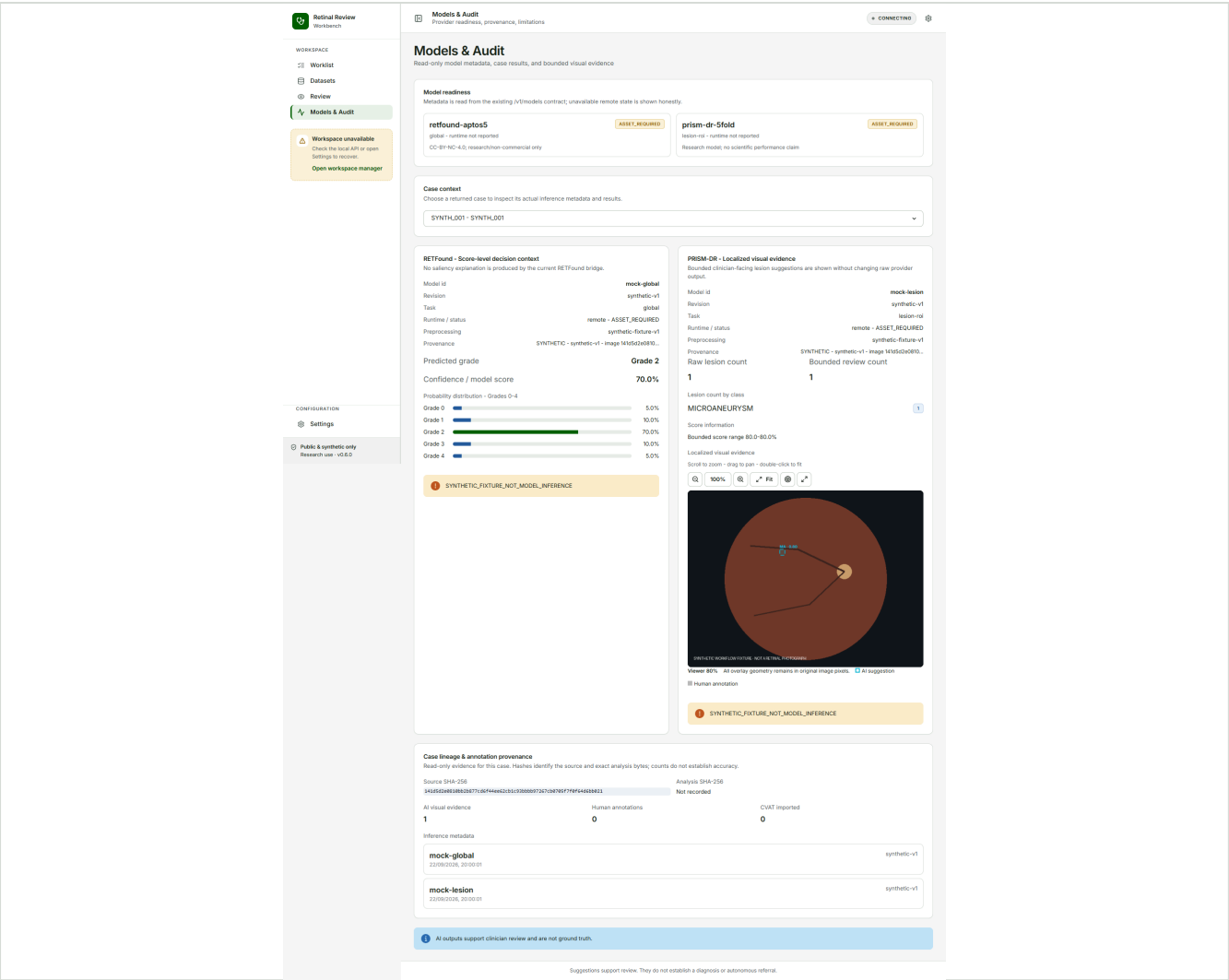

Models and audit

Safety note

Model output supports clinician review. A confidence value is a model detection score, not a calibrated probability. An empty PRISM result does not prove absence of lesions. This software is not autonomous diagnosis or referral software.